



Philips Sonicare AirFloss
Interdental - Rechargeable

- w/ 3 Nozzles
- 1 Nozzle holder
- Travel case

AirFloss



HX8255/02

Don't floss? Then AirFloss.

Removes plaque where brushing can't

For those who don't floss consistently, AirFloss is the easiest way to start cleaning between teeth. AirFloss can be used with mouthwash or water and features unique air and micro-droplet technology to remove plaque where brushing cannot.

Achieve a 30-second clean

- Cleans your entire mouth in just 30 seconds
- Point. Press. Clean!

Helps build a healthy habit

- An easy way to improve cleaning between teeth

Improves oral health

- Removes more plaque between teeth
- Helps improve gum health in just two weeks
- Helps prevent cavities between teeth

Easy to use [AirFloss]

- Air and micro-droplet technology
- Comfortable ergonomic handle
- Guidance tip ensures correct placement

PHILIPS
sonicare

Highlights


An easy way to start

Interproximal cleaning is very important to overall oral health. AirFloss is an easy way to clean deeper between teeth, helping to form a healthy habit. After three months of using AirFloss, 96% of inconsistent flossers continued to use AirFloss four or more days per week.

Removes more plaque between teeth

Removes more plaque between teeth

Air and micro-droplet tech

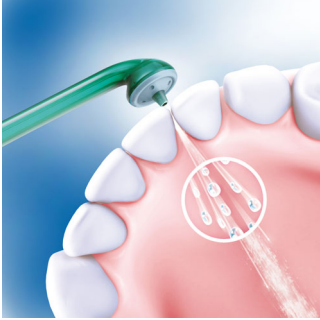

Air and micro-droplet technology

Ergonomic handle

Designed to be easy to maneuver in tight spaces, AirFloss has a comfortable ergonomic handle and a non-slip grip.

Guidance tip



The slim angled nozzle and guidance tip make finding the correct placement easy. Just slide along the tooth's surface near the gum line until the tip fits into the groove between the teeth.

Auto-burst feature



With the auto-burst feature, cleaning your entire mouth takes just 30 seconds. Simply hold down the button and move the guidance tip from space to space. It will automatically burst every second. For those who want more time between bursts, press the button once for a single burst, then move and repeat.

Easy-fill reservoir



Fill the reservoir on the handle with either mouthwash or water for a refreshing clean, then point and shoot. The reservoir holds enough liquid for two uses.

Better gum health

Sonicare AirFloss is clinically proven to improve gum health in just two weeks.

Helps prevent cavities

By gently bursting away plaque that brushing missed, Sonicare AirFloss helps prevent cavities from forming in the spaces between your teeth.

Specifications

Ease of use

- Handle: Slim ergonomic design for easy handling.
- Battery indicator: Flashing indicates low battery
- Battery Life: Up to 2 wks when fully charged
- Nozzle system: Easily snaps on and off

Technical specifications

- Battery: Rechargeable
- Operating time: Replace nozzle every 6 months
- Battery type: Lithium ION rechargeable
- Charge time: 24 hours to full charge

Power

- Voltage: Multi-voltage charger

Service

- Warranty: 2-year limited warranty

Items included

- AirFloss Handle: 1
- AirFloss Nozzle: 3
- Charger: 1 multi-voltage charger
- Nozzle Holder: 1

Cleaning performance

- Performance: Removes up to 5x more plaque*
- Health benefits: Improves gum health in only two weeks
- Auto-burst feature: Cleans entire mouth in 30 sec
- Nozzles: Best results change every 6m

Design specifications

- Color(s): White with green accents



Issue date 2013-06-18

Version: 2.1.1

© 2013 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* than manual toothbrushing alone