

iluminage
PRECISE TOUCH

USER MANUAL



IT'S THE START OF SOMETHING BEAUTIFUL.

USER MANUAL - ORIGINAL ENGLISH INSTRUCTIONS

These instructions, warnings and safety information within this user manual are intended to help ensure proper use and optimal results. It is important that you read this entire user manual before using your **PRECISE TOUCH** for the first time.

If you have questions or concerns, please contact our Customer Care department at:

www.iluminagebeauty.com

TABLE OF CONTENTS

CHAPTER 1 : INTRODUCTION

CHAPTER 2: GETTING TO KNOW YOUR **PRECISE TOUCH**

.....: 06 Device components

CHAPTER 3: WARNINGS

.....: 10 Who should not use (contraindications)

.....: 10 Warnings

CHAPTER 4: USING YOUR **PRECISE TOUCH** FOR THE FIRST TIME

.....: 14 Setting up your **PRECISE TOUCH**

.....: 15 Skin sensitivity test

.....: 15 How to use

.....: 17 Treatment Tips

.....: 18 What to expect during treatment

.....: 18 Treatment schedule

CHAPTER 5: TAKING CARE OF YOUR **PRECISE TOUCH**

.....: 22 How to clean and store your **PRECISE TOUCH**

.....: 22 Disposal/Recycling

.....: 22 Travel with your **PRECISE TOUCH**

CHAPTER 6: TROUBLESHOOTING

CHAPTER 7: WARRANTY

CHAPTER 8: MANUFACTURER INFORMATION

INTRODUCTION

CONGRATULATIONS!

PRECISE TOUCH is the only in-home compact hair removal device, using patented elōs technology, that is clinically proven to be safe and effective for all skin colours for permanent results on the body and face.

NEED HELP? Get more info at: www.iluminagebeauty.com

INTRODUCTION

Developed by the leading global aesthetic device company with the same technology used by dermatologists around the world, **PRECISE TOUCH** treatments disable the hair follicle so unwanted hair does not grow back. Gentle, compact and easy to use, get ready to experience silky, smooth skin without the hassles of shaving and waxing!




INTENDED USE OF THE PRECISE TOUCH

PRECISE TOUCH is an over the counter device intended for permanent hair reduction on the body for men and women and on women's facial hair below the cheek bone.

Skin Colour

PRECISE TOUCH is intended for use on ALL skin colours (the skin colour you were born with).

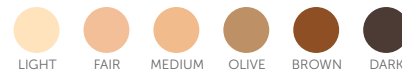

Hair Colour

PRECISE TOUCH is clinically proven to be effective on the widest range of hair colours.



Results may vary based on individual's skin colour and hair colour being treated and may require additional treatments beyond the full protocol described on page 11. For more information, please contact Customer Care at www.iluminagebeauty.com




GETTING TO KNOW YOUR **PRECISE TOUCH**


NEED HELP? Get more info at: www.iluminagebeauty.com

2

PRECISE TOUCH TOP VIEW

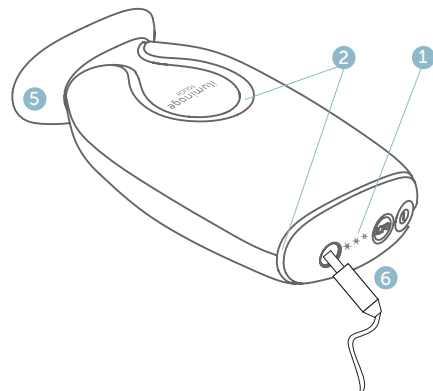

PRECISE TOUCH BOTTOM VIEW

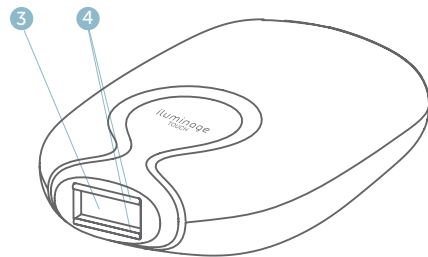

PRECISE TOUCH PACKAGE INCLUDES




PRECISE TOUCH device

Power supply

DEVICE COMPONENTS

1

CONTROL PANEL

-  Power button
-  elōs energy setting button with 3 energy level indicators: LOW- MEDIUM- HIGH
-  Power supply port

2

COOLING VENTS

Cool air is continuously circulated to keep the device cool while in use.

3

TREATMENT WINDOW

The light energy flashes out through the treatment window.

4

RADIO FREQUENCY (RF) CONTACT BARS

Both RF bars must be in complete contact with the skin for the light energy to flash.

5

PRECISE TOUCH STAND

Stand holds the device and covers the treatment window when not in use.

6

POWER SUPPLY

Connects the device to an electrical source during use.

WARNINGS

READ THE FOLLOWING IMPORTANT
SAFETY INFORMATION BEFORE USING
PRECISE TOUCH

NEED HELP? Get more info at: www.iluminagebeauty.com

WHO SHOULD NOT USE (CONTRAINDICATIONS)

DO NOT use if you have epilepsy or light-induced seizures or chronic migraines. The energy flashes might cause severe headaches and migraines.

DO NOT use if you have an active implant, such as a pacemaker, incontinence device, insulin pump, etc. Use of the **PRECISE TOUCH** has not been tested in individuals with active implants, and, therefore, the effects are unknown.

KEEP OUT OF REACH OF CHILDREN. The **PRECISE TOUCH** should not be used by anyone, or on anyone, under the age of 18. The use of the **PRECISE TOUCH** on or by children has not been studied and could result in serious damage.

The **PRECISE TOUCH** is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.

WARNINGS

DO NOT use if you are pregnant, trying to become pregnant (conceiving) or breast-feeding. The **PRECISE TOUCH** has not been tested on pregnant or breast-feeding women. Hormonal changes could increase skin sensitivity and the risk of skin damage.

DO NOT treat with the **PRECISE TOUCH** immediately after or before an active tanning session - in natural light, on a tanning

bed or when using sunless tanning products. If you are not certain whether you are tanning or sunburned it is recommended to do a skin sensitivity test (page 8) and wait 24 hours.

DO NOT use the **PRECISE TOUCH** on areas that have been sunburned as this can make the skin more susceptible to damage.

DO NOT use directly on parts of your skin where you have tattoos, dark brown or black spots (such as large freckles, birthmarks or moles). These areas may absorb too much light energy which can damage the skin.

DO NOT use on the ears, neck, nipples, genitalia or around the anus. These areas may have more sensitive skin and using the **PRECISE TOUCH** there may cause serious damage.

DO NOT use on or around the eyes, eyebrows or eyelashes. Doing so can cause permanent eye damage.

DO NOT use on a man's face, jaw or neck. Using the **PRECISE TOUCH** in these areas may result in skin damage.

DO NOT use the **PRECISE TOUCH** on any area where you may want the hair to grow back. The hair removal results may be permanent.

DO NOT use the **PRECISE TOUCH** if you are sensitive to light (photosensitive).

- Some diseases can be related to photosensitivity including porphyria, polymorphic light eruption, solar urticaria and lupus.
- Some medications can be related to photosensitivity

including: NSAIDs; steroids; antibiotics; phenothiazine; thiazide, diuretics, sulfonylureas, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), Retin-A, Accutane® and/or topical retinoids.

- Undergoing radiation or chemo-therapy treatments.

These conditions can make the skin more sensitive to light and more susceptible to damage. Consult your doctor.

DO NOT use if you have an active or known skin condition in the treatment area including acne, sores, psoriasis, eczema, sunburns, melasma, herpes simplex, open wounds or a rash. These conditions can make the skin more sensitive to light and more susceptible to damage.

DO NOT use if you have a history of skin cancer or areas of potential skin malignancies. Use of the **PRECISE TOUCH** has not been tested in areas with skin cancer, and therefore, the effects are unknown.

DO NOT hold the treatment window still over one spot on the skin and emit multiple light flashes on the skin. Doing so may generate too much heat, which can damage the skin.

DO NOT cover the cooling vents during use, as they help to keep the device cool during operation. Covering the cooling vents could result in the device overheating.

DO NOT get the **PRECISE TOUCH** damp or wet.

DO NOT use the **PRECISE TOUCH** while bathing or put it next

to a bath or sink.

DO NOT operate the **PRECISE TOUCH** with a damaged plug or cord.

For means of isolation – disconnect the device from the mains.

Use only the power supply provided with the **PRECISE TOUCH**.

DO NOT open the outer covers of the **PRECISE TOUCH** or touch any inner parts.

DO NOT attempt to repair your **PRECISE TOUCH** device.

ALWAYS keep the cord away from heated surfaces. Failure to follow these precautions can cause electrocution, electric shock or result in damage.

ALWAYS attach the AC plug from the power supply into the device first, then into the wall socket.


DO NOT use the **PRECISE TOUCH** in an environment where a spark could cause an explosion, such as near a gas source.



USING YOUR PRECISE TOUCH FOR THE FIRST TIME

NEED HELP? Get more info at: www.iluminagebeauty.com


4

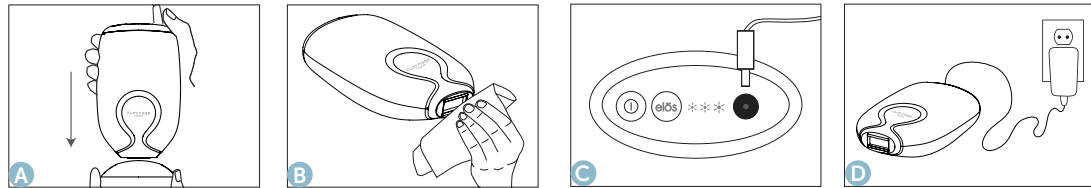

SETTING UP YOUR **PRECISE TOUCH**

- A** Lift the device from the stand. Prior to first time use, remove the protective seal sticker from the device.
- B** Clean the treatment window with a non-alcohol based wet wipe or soft tissue. This is done in order to remove any debris or dust from the treatment window.
- C** Connect the power cord into the power supply port found on the control panel of the **PRECISE TOUCH**.
- D** Connect the power cord to an AC outlet.

- E** Press the power button to turn the device on. All 3 energy level indicators will light up blue and start blinking.
- F** Press the elōs button once to select the LOW energy level. When the energy indicator lights up steady blue, the **PRECISE TOUCH** is ready for use.
- G** To increase the setting to MEDIUM or HIGH, press the elōs button until the indicator displays the desired level.
- Note:** When moving from HIGH to LOW, the device will reset. The energy indicator lights will cycle momentarily before the LOW energy indicator light appears a steady blue. Proceed treating on LOW only after the the LOW indicator light is steady.

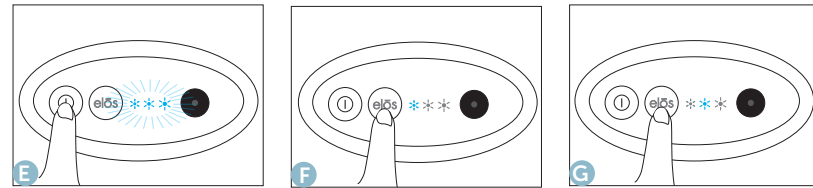

SKIN SENSITIVITY TEST

If you are using the **PRECISE TOUCH** for the first time or treating a new body area, choose a small area close to the area you intend to treat. The Skin Sensitivity Test should be performed on LOW energy level:

- Prepare your skin for the test and follow the treatment instructions.
- Wait 30 minutes and check your skin for any reaction. If the skin appears normal, then you are ready to perform a full treatment.
- If the treated spot appears red or irritated, please wait 24 hours and repeat the test. Otherwise consult your doctor.

HOW TO USE YOUR **PRECISE TOUCH**


If this is your first time using the device, be sure to conduct the Skin Sensitivity Test.

STEP 1 : PREPARATION


To prepare for your **PRECISE TOUCH** treatment:

- First remove hair from the area you want to treat.
- Clean your skin and dry it with a towel.
- This will ensure that the elōs technology is focused on the hair follicles beneath the surface for optimal efficacy.


If you have light hair, it is recommended that you wax or epilate immediately before treating with the **PRECISE TOUCH** in order to get the best results. Make sure to clean your skin of any wax or hair residue before starting treatment. Do not treat if you feel your skin is irritated from waxing or epilation.

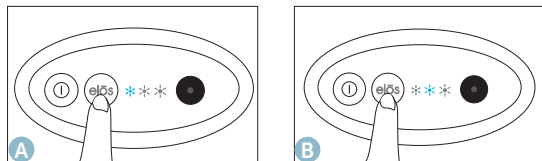

low/ faible

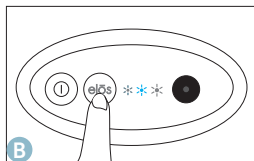

medium/ moyen

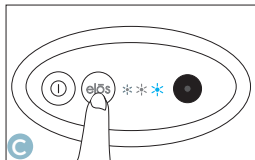


high/ élevé

STEP 2: ACTIVATE

Press the power button to turn on your **PRECISE TOUCH** device and select your energy level (LOW, MEDIUM or HIGH):

- **LOW:** Press the elōs button once to select the lowest energy level. The left energy level indicator light will turn blue (image A) and you will hear the fan inside the device.
- **MEDIUM or HIGH:** Use the elōs button to select higher energy levels for greater efficacy (image B, C). Always start treating on LOW and if you are comfortable, move to MEDIUM or HIGH energy level for greater efficacy.

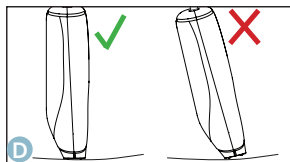

skin surface/ surface de la peau

STEP 3: TREATING YOUR BODY

Gently place the device on your skin, so that both RF contact bars are in complete contact with the skin (image D). The device will not work unless both RF contact bars are touching your skin.

- The **PRECISE TOUCH** will emit flashes of bright light - this is the elōs energy at work.
- Slowly glide the device over the treatment area, without lifting the device move to the next spot only after the treatment flash has been emitted.
- Continue gliding the device over the area so that you pass over the same area 2-3 times.



STEP 4: TREATING YOUR FACE

Prepare your treatment area (women's faces below the cheek line):

- You can shave (dark hair), wax or tweeze (light hair) before using the **PRECISE TOUCH**.
- Make sure to clean your skin of any wax or hair residue before starting treatment.
- Do not treat if you feel your skin is irritated from waxing or tweezing.
- Gently place the device on your skin, so that both RF bars are in complete contact with the skin. The device will not work unless both RF contact bars are touching your skin.
- Allow one flash of elōs energy on your first target area.
- Lift the device and place it on the the next adjacent area and wait for each flash of elōs light to be emitted before moving to the next spot.
- Repeat until you have completed 2-3 passes of the entire area.



TREATMENT TIPS

When using the **PRECISE TOUCH** on difficult to treat areas (e.g. underarms) apply the device in two directions: left/right and up/down for full coverage (image A).

- If you find the light from the elōs flash too bright, you can wear sunglasses or goggles for your comfort.
- **DO NOT** hold the treatment window still over one spot and flash multiple light flash on the skin. Doing so may generate too much heat, which can damage the skin.
- The **PRECISE TOUCH** will automatically shut off if it is not in use for more than 3 minutes. Press the power button to turn the device back on.



WHAT TO EXPECT DURING TREATMENT

During treatment you may feel a warming sensation. You should always start your treatments at the LOW energy level and if you are comfortable, increase to MEDIUM or HIGH for even greater efficacy. If your skin looks slightly red after treatment, this should disappear within 24 hours.

TREATMENT SCHEDULE

Hair removal is a gradual process. Actual results will vary from person to person and depend on the body site being treated. Hair re-growth may fluctuate over time, depending on the person, body area, hair growth cycle and whether the instructions were followed correctly.

For best results, follow the full protocol of treating the target area once a week for 7 treatments. Use your temporary hair removal solution in between treatments. After your initial protocol, use your **PRECISE TOUCH** for touch-ups or as needed.


Note: Facial areas can be treated more frequently of up to 2-3 times per week for 6 weeks if desired. Additional treatments may be needed for optimal results.



TAKING CARE OF YOUR PRECISE TOUCH

NEED HELP? Get more info at: www.iluminagebeauty.com

5

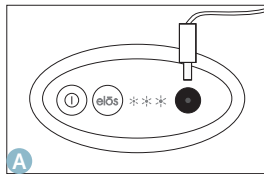


HOW TO CLEAN AND STORE YOUR PRECISE TOUCH

Turn the **PRECISE TOUCH** off by pressing the power button on the control panel. Unplug the power supply from AC outlet and power supply port (image A).

Clean the treatment window and outside surface of the **PRECISE TOUCH** using ONLY a non-alcoholic based wet wipe or soft tissue (image B).

Store the **PRECISE TOUCH** in a cool and dry place AWAY from the reach of children.

It is recommend to store your **PRECISE TOUCH** along with its stand to cover the treatment window when not in use (image C).

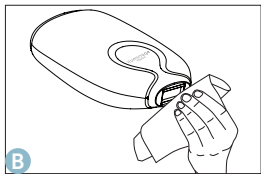

TRAVELING WITH PRECISE TOUCH

The **PRECISE TOUCH** is designed with universal voltage. That means that you can use the device all over the world. Interchangeable plugs are sold separately and vary by country. These are designed to work with the **PRECISE TOUCH** power supply. For information on purchasing, please contact us at www.iluminagebeauty.com. Make sure that the **PRECISE TOUCH** does not get wet while travelling, and do not use it in a wet environment.

DISPOSAL/RECYCLING

Please follow your local regulations for disposal of electric systems when recycling or disposing of your **PRECISE TOUCH**. If you have questions, please contact us at www.iluminagebeauty.com.



Waste Electrical and Electronic Equipment (WEEE).



TROUBLESHOOTING

NEED HELP? Get more info at: www.iluminagebeauty.com

6

EN TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
No indicator lights are on.	<ul style="list-style-type: none"> Power supply is not properly connected to the device or to an electrical source. Device wasn't turned on. 	<ul style="list-style-type: none"> Ensure the DC plug of the power supply is securely plugged into the device and the other end is connected to an AC outlet. If the power supply has been properly connected, press the power button located on control panel of the device.
All 3 energy indicator lights are blinking blue.	<ul style="list-style-type: none"> An energy level was not selected after turning on the device. 	<ul style="list-style-type: none"> Press the elōs button once to select the LOW energy level. To select MEDIUM or HIGH energy levels, use the elōs button until the indicator light displays the energy level you desire.
No energy light flashes from treatment window.	<ul style="list-style-type: none"> Device wasn't turned on. Energy level was not selected by pressing the elōs button. RF skin contact bars not in contact with the skin in treatment area. Device was left idle for more than 3 minutes and turned off. 	<ul style="list-style-type: none"> In order for the PRECISE TOUCH to emit flashes, make sure all 3 conditions (device on, energy level selected, RF skin contact bars in contact with skin) are met. If the device was idle, press the power button to reactivate the device for treatment again.
All 3 energy indicator lights are steady blue (not blinking).	<ul style="list-style-type: none"> Device flashes consumed or a technical issue. Cooling vents covered. Room temperature too high. 	<ul style="list-style-type: none"> Turn off and disconnect the device from the power supply and contact our Customer Care Specialists for assistance. Let device cool. Turn the device back on to re-attempt your treatment. Plug the cord back into the device. Move to a cooler room
3 energy indicators are blinking in a continuous cycle.	<ul style="list-style-type: none"> When the device was set to HIGH energy level the elōs button was pressed to lower the energy level to LOW. 	<ul style="list-style-type: none"> Wait for the indicator lights to stop cycling and the LOW energy level indicator is a steady blue. Then treat the desired area.

If you still have questions about your **PRECISE TOUCH** system, contact Customer Care at www.iluminagebeauty.com

WARRANTY

NEED HELP? Get more info at: www.iluminagebeauty.com

7

ILUMINAGE PRECISE TOUCH WARRANTY

NOTE: if your **PRECISE TOUCH** is under warranty with your local distributor it will prevail over the below warranty terms.

Your **PRECISE TOUCH** is warranted to be free from defects in materials and workmanship for a period of 2 years (or otherwise stated by your local distributor/retailer) from the original purchase date under normal use. This warranty extends only to the original purchaser with original proof of purchase and only when purchased from an authorised **PRECISE TOUCH** retailer or reseller. If the product should become defective within the warranty period, please contact your local retailer/distributor or go to www.iluminagebeauty.com for further assistance.

Iluminage LTD. ("Iluminage") reserves the right to replace a defective product with the most comparable product currently available. Your exclusive remedy for any defective product covered by this warranty is limited to the repair or replacement of the defective product. Iluminage's entire liability for any defective product shall in no event exceed the original purchase price of the defective product.

TO THE EXTENT PERMITTED BY LAW, Iluminage SHALL NOT BE RESPONSIBLE FOR ANY DIRECT OR INDIRECT DAMAGES RESULTING FROM THE USE OR MISUSE OF THIS PRODUCT INCLUDING ANY GENERAL, INCIDENTAL, SPECIAL, OR CONSEQUENTIAL DAMAGES - TO THE MAXIMUM EXTENT PERMITTED BY THE LAW. THE WARRANTY PROVIDED IN THIS DOCUMENT IS THE SOLE AND EXCLUSIVE WARRANTY PROVIDED BY Iluminage TO THE PRODUCT AND IT REPLACES ANY REPRESENTATION AND/OR WARRANTY INCLUDING FITNESS FOR ANY PARTICULAR USE AND MERCHANTABILITY.

This warranty does not cover products determined by Iluminage to be damaged by any of the following: accident, misuse, abuse or alteration, servicing by unauthorised personnel, use with unauthorised accessories or use other than as instructed, connecting to an incorrect current or voltage and/or any other conditions beyond Iluminage's control.



MANUFACTURER INFORMATION

NEED HELP? Get more info at: www.iluminagebeauty.com

8


LABELS AND SYMBOLS

The **PRECISE TOUCH** device shows information pertaining to the device's part number, serial number, manufacturer information and the electrical and regulatory standards the **PRECISE TOUCH** complies with.

This device is not suitable for use in the presence of a flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.



Illuminage Ltd.
Yokneam Illit 2069200, IL

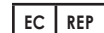

mdi Europa GmbH
Langenhagener Str. 71
D-30855 Hannover-Langenhagen

DEFINITION OF SYMBOLS



WEEE - Waste Electrical and Electronic Equipment



Follow operating instructions



Class II equipment



Do Not use in wet environment



Degree of protection against electric type BF applied part



CE Mark

PRECISE TOUCH SPECIFICATIONS*

PRECISE TOUCH

Light spectrum	550-1200 nm
Max optical energy	up to 4 J/cm ²
RF wave frequency	6.78 MHz
Max radio frequency power	up to 5W avg.
Flash light repetition rate	up to 1 Hz (1 flash every 1 second)
Electrical rating	100-240Vac, 0.9Amax, 50-60Hz
Dimensions	4 x 9 x 15 cm (W x L x H)
Weight	0.25 kg

*Patented technology US67028208 - US5405368

APPLICABLE STANDARDS

SAFETY STANDARDS

EN 60335:1 2002 (A1:2004; A2:2009)	IEC 60601-2-57
EN 60335-2-27: 2009	ETL CLASSIFIED
IEC 62471:2006	CONFORMS TO UL STD.
IEC 60601-1, 3rd edition	UL 60601-1 CERT TO CSA STD.
IEC 60601-1-11	C22.2 No. 601.1

EMC STANDARDS

EN 61000-6-1:2007
EN 61000-6-3:2007
IEC 60601-1-2:2007

ACCEPTABLE ENVIRONMENTAL CONDITIONS

STORAGE/TRANSPORT

Temperature	: 14°F - 131°F/ -10°C - 55°C
Humidity	: Up to 90% at 131°F
Pressure	: 700-1060 hPa

OPERATING

Temperature	: 50°F-86°F/ 10°C-30°C
Relative Humidity	: Up to 80% at 99°F
Altitude	: Up to 13,000ft

